
REPORT OF ANALYSIS - PHYSICS, PARTICLE & SURFACE ANALYSIS

Sample: MRa389A
Order title/Product: Graphene / Raman microscopy
Project: Graphene
Account: Curia Germany GmbH
Customer: Dr. Thorsten Busch
Order no.: 25-049538-4
Remarks: Raman measurement on request of the customer

Method:

Raman spectroscopy/microscopy using Senterra I Raman microscope (Bruker).

Procedure:

Objective: Raman spectroscopy
Field of Application: Powder
Principle: The substance is applied to a microscope slide and examined under the microscope. The objective is selected and focused to examine an exact location using Raman. At this position, several points are examined with the laser and a corresponding number of Raman spectra are recorded.
Reference:

Experimental Details and Result

Equipment: Senterra I Raman microscope (Bruker)
Magnification: 20 x
Number of spectra per position: 8
Evaluation: Averaging over several spectra

Result:

Absorption band	2D band	G band	D band	D/G ratio
Integral	n.d.	216,4	236,7	1,094

Spectra see attachment.

Electronic document valid without signature!

25-049538-4: MRa389A 20x

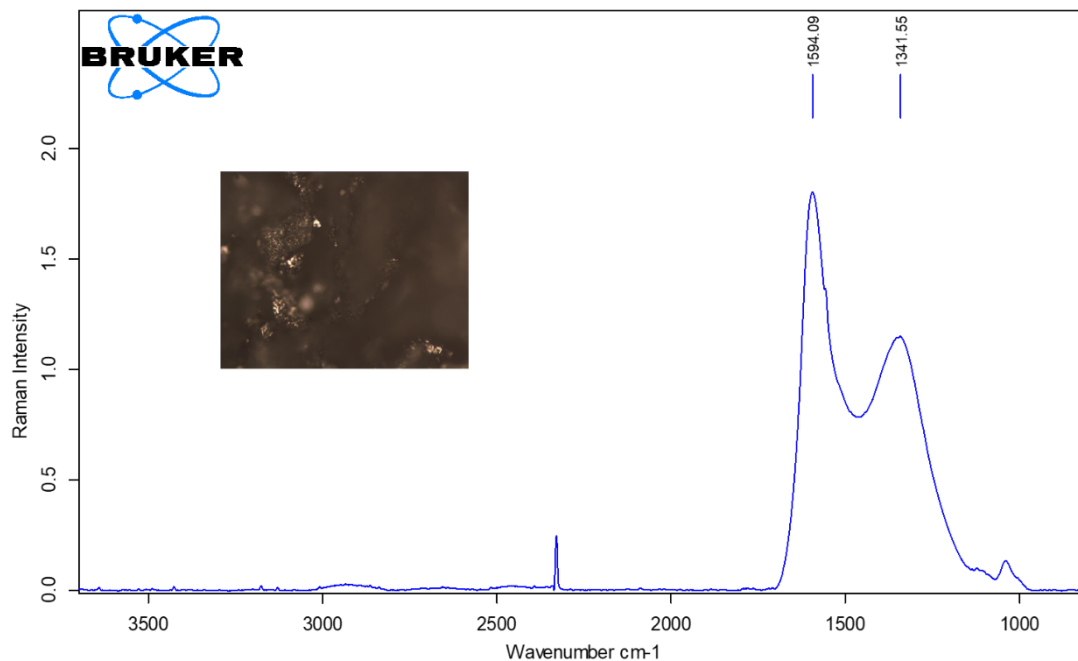


Figure 1: MRa389A, blue curve: average spectrum over 8 recordings

Electronic document valid without signature!